

### Incubating Rocking & Waving Shakers



## Save Incubator Space with Unique Benchtop Incubating Rocking and Waving Shakers

OHAUS Incubating Rocking and Waving Shakers are innovatively designed to incubate samples providing accurate and repeatable results. Rocking shakers provide a controlled see-saw like motion, while waving shakers offer a smooth, low-foaming three-dimensional wave motion. All models feature independent displays with intuitive touchpad controls. Both models include safety features that protect both the user and samples.

#### **Standard Features:**

- Electronic control allows simple tilt angle, speed and temperature adjustments while the unit is still running to achieve optimal mixing.
- Achieve reproducible results you can trust. Precisely adjust the tilt angle in 1° increments, vary the speed by 1 rpm and set temperature in 1°C increments with intuitive touchpad controls.
- Innovative speed ramping feature slowly increases speed to desired setpoint to avoid splashing. Smooth rocking or waving motion is maintained with unique, state-of-the-art motor and electronics.

# Incubating Rocking & Waving Shakers

Model	ISRK04HDG	ISWV02HDG	
Shaker Type	Rocking Shaker	Waving Shaker	
Control	Digital		
Motion / Tilt Angle	Rocking, 0 – 15°	Waving, 0 – 20°	
Speed Range	1 – 50 rpm	1 – 30 rpm	
Speed Accuracy	± 1 rpm		
Temperature Range	Ambient +5 °C − 65 °C		
Temperature Uniformity	± 0.5 °C at 37 °C		
Timer	1 second – 160 hours		
Capacity	4.5 kg	2.3 kg	
Working Environment	5 °C – 40 °C, 80% Relative Humidity, non-condensing		
Drive System	Stepper Motor		
Included Tray (L × W)	ded Tray (L $\times$ W) 254 mm $\times$ 191 mm 235 mm $\times$ 184 n		
Tray Construction	Aluminum		
Overall Dimensions (L $\times$ W $\times$ H)	419 × 257 × 274 mm		
Inner Dimensions (L $\times$ W $\times$ H)	273 mm × 197 mm × 970 mm	273 mm × 197 mm × 860 mm	
Net Weight	5.7 kg	6.2 kg	
Power	230V, 1.87A, 50/60Hz		
Power Consumption	420 W		

#### **Other Standard Features and Equipment**

Spill-resistant design, caution hot indicator, detachable 3-wire cord and plug (included)

#### Compliance

- Product Safety: EN 61010-1, EN 61010-2-010, EN 61010-2-051
- Electromagnetic Compatibility: EN 61326-1 Class A, Industrial Environments
- Compliance Marks: CE; RCM; TÜV SÜD


#### **Applications**

Hybridizations, cell cultures, blotting techniques, staining and destaining gels

#### **Accessories**

I	Incubating Rocking and Waving Shaker Accessories				
	Item No.	Description	ISRK04HDG	ISWV02HDG	
	30400141	Dimpled Mat, Incubating Rocker ISRK04HDG	•		
Г	30400143	Dimpled Mat, Incubating Waver ISWV02HDG		•	

#### **Dimensions**

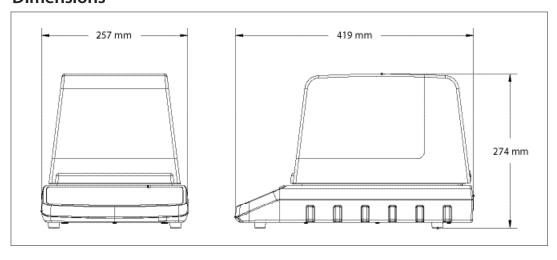

### OHAUS Asia/Pacific Headquarters

7F, Block 33, 680 Guiping Road Shanghai 200233 China

e-mail: APmarketing@ohaus.com Tel: +86 21 64855408

#### www.ohaus.com

The management system governing the manufacture of this product is ISO 9001:2015 certified.

SAP Literature #



80775376\_A 20190503 © Copyright OHAUS Corporation