

CIRCULATED WATER BATHS

(AMBIENT +5°C TO +100°C)

Image: TWBC-24 shown with optional TU3



Thermoline Digital Circulated Water Baths

- 12, 24 and 48 litre models available
- Bath constructed of high grade stainless steel
- Insulated walls to reduce heat loss
- Seamless stainless steel tank with radial corners
- Carry handles for easy moving
- Custom made baths available on request
- Australian made

Instrument Choice

Phone: 1300 737 871

Fax: 1300 882 503

Email: customer-service@instrumentchoice.com.au

www.instrumentchoice.com.au

 **Thermoline**
S C I E N T I F I C

HIGHER STANDARDS IN SCIENCE

PRODUCT DETAILS

Thermoline's quality range of stainless steel laboratory water baths available in 12, 24 and 48 litres capacity. All baths are constructed of stainless steel, and are insulated to reduce heat loss. Accurate temperature control is available when used together with the TU series heater circulators.

Optional Accessories

Polycarbonate gable lids to suit 12 and 24 litre models

Stainless steel gable lid to suit 48 litre model

Stainless steel flat lids available to suit all models


TU Series heater circulators to allow heating from ambient +5°C to 100°C

- For more information on Thermoline's range of Unistat Heater Circulators, please refer to our "Digital Heater Circulators" page.
- Information on Thermoline's Uncirculated Baths can be found on our "Digital Uncirculated Baths" page.

SPECIFICATIONS

Model	TWBC-12	TWBC-24	TWBC-48
Capacity	12 Litres	24 Litres	48 Litres
Internal Dimensions (WxDxH mm)	160x230x200	370x300x200	490x500x200
External Dimensions (WxDxH mm)	330x270x230	530x330x230	650x530x230
Optional Gable Lids	PCLC-12 (polycarb)	PCLC-24 (polycarb)	SSGLC-48 (stainless)
Optional Flat Lids	SSFLC-12 (stainless)	SSFLC-24 (stainless)	SSFLC-48 (stainless)
Weight (kg)	6	8	15

Image: TWBC-48 shown with optional TU1

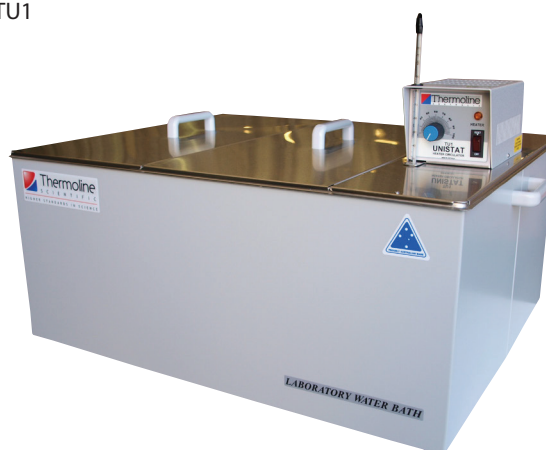

Thermoline Scientific have been manufacturing and distributing high quality laboratory and scientific testing equipment since 1970. Over this time, Thermoline has grown to be a leading brand in the science industry, with our equipment being used in small and large Hospitals, Universities and Research Laboratories across Australia and the Asia Pacific region.

We're proud to say that we are 100% Australian owned and operated.

Instrument Choice

Phone: 1300 737 871

Fax: 1300 882 503

Email: customer-service@instrumentchoice.com.au

www.instrumentchoice.com.au

CIRCULATED WATER BATHS

(AMBIENT +5°C TO +100°C)

Image: LWB-311DS with Gable Lid



Labtech Digital Circulated Water Baths

- 11 and 22 litre models available
- Unique built in water circulation pump
- Bath constructed of high grade stainless steel
- Economic and easy to operate
- Seamless stainless steel tank with radial corners
- Flat and gable lids available for both models

Instrument Choice

Phone: 1300 737 871

Fax: 1300 882 503

Email: customer-service@instrumentchoice.com.au

www.instrumentchoice.com.au

 **Thermoline**
S C I E N T I F I C
H I G H E R S T A N D A R D S I N S C I E N C E

PRODUCT DETAILS

These two circulated water baths are powerful and economic. Both models feature a unique in-built circulation water pump and are fitted with digital PID control systems which provide accurate temperature control. A front touch membrane keypad allows users to make easy temperature, timer and alarm adjustments. Temperature auto-tuning function is supported for user calibration.

All baths are supplied standard with a flat stainless steel lid with the option of a clear gable lid.

Optional Accessories

Gable lids available for all models

SPECIFICATIONS

Model	LWB-311DS	LWB-322DS
Capacity	11	22
Bath Dimensions (WxDxH mm)	300x240x150	500x295x150
Overall Dimensions (WxDxH mm)	327x265x280	537x327x280
Temperature Range	Ambient +5 °C to 100 °C	
Controller	Digital PID Control System	
Safety	Over temperature protection with audible alarms	
Pump	9 Litres / Min	
Timer	99 Hr 59 min	

Thermoline Scientific have been manufacturing and distributing high quality laboratory and scientific testing equipment since 1970. Over this time, Thermoline has grown to be a leading brand in the science industry, with our equipment being used in small and large Hospitals, Universities and Research Laboratories across Australia and the Asia Pacific region.

We're proud to say that we are 100% Australian owned and operated.

Instrument Choice

Phone: 1300 737 871

Fax: 1300 882 503

Email: customer-service@instrumentchoice.com.au

www.instrumentchoice.com.au